**Over-the-Wire LeMaitre® Valvulotome**

English – Instructions for Use

**Over-the-Wire LeMaitre® Valvulotom**

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**Valvulotome Over-the-Wire LeMaitre®**

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**Valvulotomo over-the-wire (OTW) LeMaitre®**

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**Valvulotomo guiado LeMaitre®**

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**Valvulótomo através de um fio guia LeMaitre®**

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**STERILE EO**

## Description

The Over-the-Wire LeMaitre® Valvulotome is used to cut saphenous vein valves during bypass procedures. The Centering Hoops keep the head of the Over-the-Wire LeMaitre Valvulotome centered in the vessel and prevents the valve-cutting blades from damaging the vessel wall. The size of the Centering Hoops and Cutting Blades adjust to the internal diameter of the vein as the Over-the-Wire LeMaitre Valvulotome is being drawn through the vessel. The LeMills Valvulotome is provided as an accessory device to cut the last valve set as needed.

Kit includes LeMills Valvulotome.

## Indication for Use

The Over-the-Wire LeMaitre Valvulotome is indicated for cutting saphenous vein valves during in-situ bypass.

The LeMills Valvulotome is intended to cut venous valves.

## Contraindications

1. pregnant women
2. patients with allergy to any of the device materials
3. endarterectomy procedures
4. thrombolysis procedures
5. vein stripping procedures
6. embolectomy procedures, or
7. vessel dilation procedures

## Warnings

1. Do not use if the package or device is damaged.
2. Do not use if packaging is opened outside a sterile environment.
3. The Over-the-Wire Valvulotome and the LeMills Valvulotome are to be used by a qualified physician.
4. Do not insert the Over-the-Wire Valvulotome into a vessel or extract from a vessel in the open position.
5. The Over-the-Wire Valvulotome must be used with an 0.035" Guidewire.
6. Do not open or close the Over-the-Wire Valvulotome while in a coiled configuration.
7. The Over-the-Wire Valvulotome must be flushed with saline or heparinized saline only.
8. Do not pass the Over-the-Wire Valvulotome or LeMills Valvulotome through a vessel that has undergone synthetic grafting or contains implants.
9. Do not rotate the Over-the-Wire Valvulotome in a vessel.
10. A maximum of 8 passes with the Over-the-Wire Valvulotome may be completed.

## Precautions

1. United States Federal and other law restricts this device to sale on or by the order of a physician.
2. Do not use device past the expiration date printed on the labeling.
3. These devices are single use only. Do not reuse, reprocess, or re-sterilize (reference "Re-sterilization/Repackaging" section of this document for further details).
4. Keep any dissected portion of the vein straight and avoid twisting the vein during valvulotomy.
5. Use caution when advancing the guidewire past the proximal anastomosis.
6. Use caution when tracking the Over-the-Wire Valvulotome in thrombophlebitic veins.
7. Do not advance the Over-the-Wire Valvulotome with the blades in the open position.

## Potential Complications


- Vessel damage
- Vessel occlusion/ stenosis
- Hematoma
- Hemorrhage
- Thrombus formation
- Wound infections
- Erythema
- Myocardial infarction
- Amputation
- Death

## Additional Required Items


1. .035" Guidewire
2. Radiopaque Tape

## Specifications

Over-the-Wire Valvulotome:

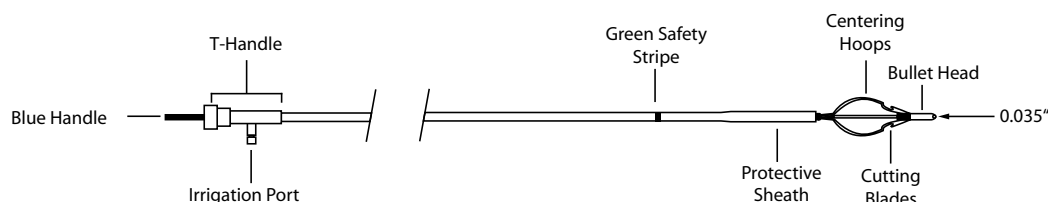LeMills Valvulotome:

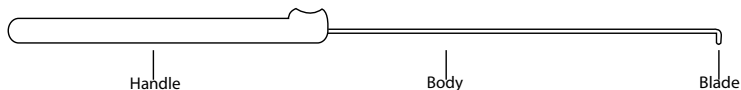Catalog Number	Useable Length	Closed Catheter Outer Diameter	Maximum Open Cutting Blade Diameter	Maximum Open Centering Hoop Diameter	LeMills Valvulotome
1005-00, e1005-00	100 cm	2.0 mm	6.0 mm	9.5 mm	Included
1006-00, e1006-00	40 cm	2.0 mm	6.0 mm	9.5 mm	Not Included

**To Open Package**

1. Open box and remove sealed sterile unit.
2. Open sealed sterile unit and present contents to personnel in sterile environment.
3. Carefully remove lids, as applicable.  
**Caution: Remove product from tray with care to ensure centering hoops are not damaged.**
4. Place in basin filled with sterile, heparinized saline.

**Pre-Use Check**

Inspect centering hoops and cutting blades for damage and alignment prior to use.

**Instructions for Use**

1. Perform an intra-operative venogram using Radiopaque Tape; mark the fistulae.
2. Excise the Sapheno-Femoral Complex.
3. Excise first valve set under direct vision using standard technique.
4. Perform the proximal anastomosis.
5. Ensure that the leg is fully extended and pass the 0.035" guidewire well beyond the treatment area.
6. Flush the Over-the-Wire Valvulotome (in the open position) by removing Irrigation Port cap, blocking guidewire lumen on the Blue Handle, and flushing with sterile heparinized saline. Replace cap upon completion.
7. Sheath the Centering Hoops within the Protective Sheath by advancing the T-Handle until Cutting Blades are fully enclosed and a positive stop is detected.
8. Load the device over the guidewire carefully.
9. Advance into the vein to the treatment area; 2 cm to 3 cm distal to the anastomosis.
10. Expose the Cutting Blades and Centering Hoops by maintaining the position of the Blue Handle and retracting the T-Handle until a positive stop is detected.  
**Note: Failure to maintain the position of the Blue Handle while exposing cutting blades may cause damage to the anastomosis.**
11. Slowly retract the Over-the-Wire Valvulotome in the open position to disrupt valve sets.  
**Warning: Do not rotate the Over-the-Wire Valvulotome.**
12. When the Green Safety Stripe is visible, sheath the Cutting Blades and Centering Hoops into the Protective Sheath by maintaining the position of the Blue Handle and advancing the T-Handle.  
**Warning: Failure to sheath Cutting Blades and Centering Hoops prior to retracting the Over-the-Wire Valvulotome from the vessel may cause damage.**  
**Note: If the Over-the-Wire Valvulotome is set aside for further use within the same procedure, device must be stored in the open position in a basin of heparinized saline. Flush device as noted in step 6 prior to reintroduction.**
13. Ligate the venous tributaries previously located and marked.  
**Note: Do not ligate the last major tributary until after the distal anastomosis is complete in order to allow blood run off.**
14. Pass the device through the vessel an additional 1-2 times after ligating the venous tributaries.
15. Steps 9 through 12 may be repeated as necessary. Use the LeMills Valvulotome to excise last valve set following standard technique as needed.  
**Warning: A maximum of 8 passes with the Over-the-Wire Valvulotome may be completed.**
16. Confirm free flow through saphenous vein using doppler.
17. Complete distal anastomosis to the appropriate arterial segment.
18. Perform a completion arteriogram using radiopaque tape to locate defects.

**Storage**

Store in a cool, dry place.

**Re-sterilization/Repackaging**

These devices are single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only. If for any reason this device must be returned to LeMaitre Vascular, place it in its original packaging and return it to the address listed on the box.

**Safe Handling and Disposal**

This device is single-use and disposable device. Do not implant. Please return the used device only at the time that the device has not performed as intended or the device is related to an adverse event. In other situations, the device should not be returned but disposed accorded to local regulations.

If serious medical incidents should arise during use of this medical device, users should notify both LeMaitre Vascular and the Competent Authority of the country where the user is located.

This product contains no heavy metals or radioisotopes, and is not infectious or pathogenic. Please consult local regulations to verify proper disposal of sharps.

**Cleaning:**

1. Devices considered necessary to return should be cleaned using one of the following:
  - a) Sodium hypochlorite solution (500-600 mg/l), or
  - b) Peracetic acid solution with subsequent ultrasonic treatment

2. Devices should then be decontaminated with either:
  - a) 70% solutions of ethanol or isopropanol for a minimum of 3 hours or,
  - b) Ethylene oxide gas
3. Devices should be completely dried prior to packaging.

**Packaging:**

1. Cleaned devices should be sealed and packed in a manner that minimizes potential for breakage, contamination of the environment or exposure to those handling such packages during transit. For devices capable of penetrating or cutting skin or packaging material, the primary packaging must be capable of maintaining the product without puncture of the packaging under normal conditions of transport.
2. The sealed primary container should be placed inside watertight secondary packaging. The secondary packaging should be labelled with an itemized list of the contents of the primary receptacle. Cleaning methods should be detailed if possible.
3. Both primary and secondary packaging of cleaned, decontaminated single-use disposable devices should be labelled with an ISO 7000-0659 Biohazard symbol.
4. Primary and secondary packaging must then be packaged inside an outer package, which must be a rigid, fiberboard box. The outer shipping container must be provided with sufficient cushioning material to prevent movement between the secondary and outer containers.
5. Shipping paper and content marking for the outer shipping container are not required.
6. Packages prepared in the above manner may be shipped to:

LeMaitre Vascular  
Attn: Complaint Lab  
63 Second Avenue  
Burlington, MA 01803, USA

**Limited Product Warranty; Limitation of Remedies**

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device and that this device is suitable for the indication(s) expressly specified in these instructions for use. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERMS INCLUDE LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

These limitations do not apply to consumers in Australia or to the extent they are precluded by local law in any other jurisdiction.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre Vascular GmbH