## **Periscope™ Dissector**

Instructions for Use - English

## **Periscope™ Dissektor**

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**STERILE EO**

### **Device Description**

The Periscope™ Dissector consists of two principal components, a shaft with an angled spatula at the distal end and a hollow flexible needle running the full length of the shaft to the spatula. The shaft and angled spatula are constructed of stainless steel while the hollow needle is made of nickel titanium. The inner lumen of the shaft extends from the handle to the palm of the spatula.

The Periscope Dissector is available in one size. The useable length of the spatula/shaft is approximately 40 cm and the hollow needle can accommodate up to a 0.018" diameter guidewire.

### **Indications and Usage**

The Periscope Dissector is intended for remote intraluminal stripping of peripheral blood vessels during vascular reconstruction and to provide a conduit for guidewire placement.

General indications for peripheral vessel reconstruction include disabling claudication, ischemic pain at rest, ischemic skin lesions or gangrene.

Indications for endarterectomy specifically as a procedure for vessel reconstruction include localized lesions and heavy circumferential plaque deposits.

### **Contraindications**

- Contraindications include, but are not limited to:
- Any evidence of aneurysmal change or severe calcification may contraindicate endarterectomy since continued aneurysmal degeneration of the adventitia may be expected.
- In some cases, the vessel wall layers may be quite coherent and the plaque segment may be difficult to separate from the vessel wall. When the vessel is in this state of disease, it is often advisable to abandon endarterectomy and consider an alternative method of vascular reconstruction.
- Use of this device is contraindicated in the coronary, internal carotid arteries and other vessels of the central circulatory system.

### **Warnings and Precautions**

The Periscope Dissector may only be used by physicians trained in vascular surgery, specifically endarterectomy procedures. The techniques described herein are offered to clinicians as a guide.

#### **Periscope Dissector - Handling Precautions**

Always inspect the sterile packaging before use. If the packaging is damaged, do not use the product and contact your LeMaitre Vascular representative.

Do not re-sterilize the Periscope Dissector or the hollow needle.

The Periscope Dissector must be used prior to the "Use By" date printed on the package.

#### **Periscope Dissector - Use Precautions**

Always use fluoroscopy when advancing or manipulating endovascular devices. Do not attempt to forcibly advance the Periscope Dissector.

Confirm that the Periscope Dissector is properly oriented and the tip of the hollow needle remains retracted within the shaft while advancing the device into the vessel.

Use of the Periscope Dissector may lead to vessel perforation requiring additional interventional procedures (e.g., bypass graft surgery of the artery).

### **Potential Adverse Events**

Adverse events potentially associated with the use of the Periscope Dissector (e.g., all complications associated with PTA, remote endarterectomy or atherectomy procedures,) include but are not limited to the following:

- Incisional/puncture site inflammation, hematoma or bleeding

- Wound infection (incisional)
- Vessel injury, dissection, perforation, rupture or pseudoaneurysm
- Vessel spasm
- Emboli, thrombosis or acute occlusion
- Fever, drug reactions or allergic reactions to contrast medium
- Conversion to bypass
- Limb loss due to distal embolization
- Renal failure
- Death

### **How Supplied**

**STERILE:** Each Periscope Dissector is sterilized using Ethylene Oxide (EO) and is supplied sterile for single use only. Do not reuse or resterilize. Resterilization may adversely affect proper mechanical function and could result in patient injury.

#### **CONTENTS:**

One (1) Sterile Periscope Dissector with an Integrated Hollow Needle

One (1) Instructions for Use

**STORAGE:** Store in a cool, dry place.

### **Inspection**

Carefully inspect the sterile package for damage or defects before opening. Do not use product after the "Use By" date on the package. If the integrity of the sterile package has been compromised or the package or product is defective, do not use the product and contact LeMaitre Vascular for return information.

### **Recommended Material**

Prepare the following material using sterile technique:

- 10cc syringe filled with sterile saline
- Up to a 0.018" guidewire of appropriate length
- 10cc syringe filled with sterile contrast solution

### **Directions For Use**

Specific patient preparation and procedural techniques are determined by the physician. One suggested technique is as follows:

1. Patient preparation and administration of anticoagulation/antiplatelet regimen is left to the discretion of the physician.
2. Use sterile surgical precautions and standards.
3. Using appropriate anesthesia, gain access at the appropriate site through an arteriotomy to expose the atheroma.
4. Make a slanting, bevel-type incision of the proximal atheroma to start a cleavage plane between the atheroma and surrounding media.


*Note:* *The following procedure steps should be performed under fluoroscopic guidance.*

5. Insert the Periscope Dissector into the cleavage plane of the vessel. Using a gentle pushing motion around the atheroma core, advance the device to a point distal to the occluded segment.
6. Under fluoroscopic guidance, position the Periscope Dissector at the end of the dissection plane.
7. While advancing the needle, flush the needle with contrast medium under fluoroscopy to verify proper penetration of the intimal layer into the true lumen of the vessel.
8. Insert a 0.018" guidewire into the Periscope Dissector needle and advance until the distal tip of the wire is sufficiently anchored into the true vessel lumen. Anchoring of the guidewire will allow retraction of the Periscope Dissector while maintaining guidewire position.

*Note:*

1. *Do not use polymer covered hydrophilic guidewires.*
2. *Optimum guidewire performance is achieved with the use of flex wires or wires with very long flexible tips.*

9. Retract the needle and remove the Periscope Dissector from the patient while maintaining guidewire position.
10. Following completion of the remote endarterectomy procedure, maintain guidewire position to facilitate additional interventional procedures as required.
11. Close arteriotomy using standard procedure.
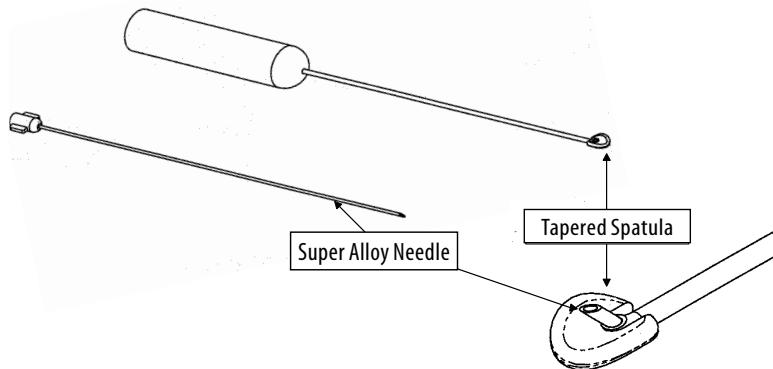
#### **Resterilization/Re-use**

This device is single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only.

#### **Limited Product Warranty; Limitation of Remedies**

LeMaitre Vascular, Inc. warrants that reasonable care has been used in the manufacture of this device. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERM INCLUDES LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. LeMaitre Vascular makes no representation regarding the suitability for any particular treatment in which this device is used, which determination is the sole responsibility of the purchaser. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY, OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

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LeMaitre is a registered trademark, and Periscope is a trademark, of LeMaitre Vascular, Inc.
©2018 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH