AlboGraft® Polyester Vascular Graft

Instructions for Use – English

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Instructions for Use – English

Description

AlboGraft® Vascular Grafts are made of synthetic material and designed to replace sections of damaged or malfunctioning arteries. They are made of polyester (polyethylene terephthalate, PET) thread woven into a seamless tube. In response to a range of surgical indications, AlboGraft Vascular Grafts are offered in two designs: double velour knitted fabric and double velour woven fabric. The knitted grafts are designed with a run-proof structure to reduce the risk of fraying or wearing down at their ends. The velour grafts have low-profile loops on their endoluminal surface to avoid any lumen reduction, and high profile loops on their outer surface to promote graft anchoring into the surrounding tissues. All AlboGraft grafts are crimped in parallel rings so that their tubular shape is maintained without kinking.

AlboGraft Vascular Grafts are available with removable external spiral reinforcement (ASC models) made of a radiopaque biocompatible polypropylene thread, allowing for easy identification of the prosthesis with x-ray. The external spiral reinforcement is removable, facilitating the creation of anastomoses to the vessel.

Note: AlboGraft Vascular Grafts with Removable External Spiral (ASC Models) are not available for sale in the United States or Canada.

Guidelines on the outer surface of the graft facilitates orientation during implantation.

The AlboGraft Vascular Grafts impregnated with collagen reduce the permeability to mitigate the problem of hemorrhaging during implantation so that no preclotting is necessary. The process of impregnation using bovine collagen maintains both the original structure of the material, and the structural characteristics of the prosthesis, i.e. flexibility and softness. The collagen is cross-linked by formaldehyde steaming to ensure a balanced reabsorption rate by the host tissues.

Available Models

AlboGraft Vascular Grafts are available in a wide range of models, types and sizes. The graft can be classified into models according to fabric characteristics, each model being of one or more types (straight "I" or bifurcated "T") and sizes (diameter, length). Each graft is identifiable by a model number and a lot/serial number. The 3-letter code identifies the fabric and presence or absence of collagen (see table below); the lot and serial number serves to retrieve information from the files of the LeMaitre Vascular Quality Assurance System on the graft's manufacturing and quality control data. To facilitate the tracking of the product by the operator all product identification data is given on self-adhesive labels in each product package, which can be transferred to the patient's clinical record.

Type	Code	Color Code
Double Velour Woven Fabric Impregnated With Collagen	ATC	light blue
Double Velour Knitted Fabric Impregnated With Collagen	AMC	gold
Double Velour Knitted Fabric Impregnated With Collagen and Spiral Supported*	ASC*	Blue

Consult the LeMaitre Vascular catalog for a detailed listing of models and types available.


Technical Characteristics of the AlboGraft Prostheses

Diameter

The diameter is defined as the value of the nominal internal diameter of the undilated prosthesis. The value is obtained with a conical gauge (anulometer) having a 10° angle from the apex. The prosthesis is gently slipped over the instrument without applying any pressure (see diagram).

D120 = Diameter of the prosthesis under pressure

** Not available for sale in the United States or Canada.*

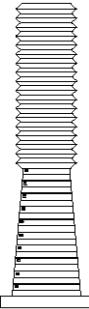This value indicates the internal diameter of the prosthesis when exposed to an internal pressure of 120 mmHg; it equals the average natural pressure on the prosthesis when unclamped. The diameter is higher than the nominal value as a result of two different concomitant factors:

- distention of the pleated folds;
- radial dilation caused by blood flow pressure (its effect on a knitted prosthesis is more sensible due to the high flexibility and elasticity of the device).

AlboGraft ASC Models*: The external spiral support ensures the radial dilation of the prosthesis remains negligible.

Water Permeability

Permeability is a chief property of the polyester fabric and depends on the porosity of the material and impregnation process. The permeability is expressed as the amount of water that can pass through 1 cm² of the prosthesis at 120 mmHg over one minute. The table below compares the maximum claimed permeability (ml x cm⁻² x min⁻¹) for each graft type.

Type	Double Velour Knitted Fabric with Collagen	Double Velour Woven Fabric with Collagen
ml x cm ⁻² min ⁻¹	0.1	0.1
3-digit code	AMC/ASC	ATC

Usable Length

This measurement indicates the minimum length of the prosthesis under tension during implantation. The minimum useable length is determined under tension with the applied weight for each type (see table below).

Type	Diameter	Weight (g)
Double Velour Knitted Fabric With Collagen (ASC)*	ALL	125
Double Velour Knitted Fabric with Collagen (AMC)	ALL	500
Double Velour Woven Fabric with Collagen (ATC)	<24 mm	500
Double Velour Woven Fabric with Collagen (ATC)	≥24 mm	750

Packaging

AlboGraft Vascular Grafts are supplied sterile through a double sterile barrier. The internal container is externally sterile provided that the external container has been neither damaged nor opened. The external container is only internally sterile: it must not be brought into a sterile environment.

Warning: The collagen-impregnated AlboGraft Vascular Graft (AMC, ATC, and ASC*) double container is provided in a vacuum sealed foil pouch. This pouch protects the graft during storage and should be opened only immediately prior to use.

The labeling is color coded to identify the fabric characteristics and collagen impregnation of the model contained inside (see table within "Available Models").

Indications

The AlboGraft Knitted and Woven Polyester Grafts are indicated for use in the replacement or repair of arteries affected with aneurismal or occlusive disease.

Instruction For Use

Choice and Preparation of the Graft

Choose the graft type and size best suited to the vascular section to be replaced or supplied through a bypass.

* Not available for sale in the United States or Canada.

Carefully examine the double packaging enclosing the selected graft to verify its integrity, the absence of holes, cuts or perforations, and the expiration date. The sterility of the graft is guaranteed only if the container is perfectly intact. Make sure that the period of sterility has not expired, otherwise the sterility of the graft is not guaranteed. Make sure that the information on the labels relative to the model, type, size, and lot number is correct. If it is not, do not use the graft. Remove Tyvek lid from external container and introduce inner container into sterile environment using aseptic procedures.

CAUTION: do not bring the external container into the sterile environment.

Remove the graft using an aseptic procedure.

Preparation Vascular Grafts

The collagen impregnated vascular grafts can be used as supplied or after immersing it in a sterile saline solution just before implantation, to make them flexible and easy to handle. During the preparation and the operation, care should be taken not to crush or excessively expand the graft. Correct handling helps avoid the risk of damaging the collagen impregnation on the graft walls.

AlboGraft Removable Spiral Reinforcement Prosthesis (ASC Models)*

Cut the graft at the desired length. Remove the spiral support at the ends and cut close to the graft surface with sharp surgical scissors in order to obtain the desired spiral free length to create the anastomosis, thus preventing the spiral support from interfering with suturing. Avoid clamping the graft on reinforced area.

Surgical Technique

Considering how complex a surgical implantation of vascular grafts is and how many factors can influence the patient's condition, it is left to the surgeon to define suture technique and graft type as well as the therapy to adopt before, during and after the operation.

The following precautions should; however, be considered:

- for suturing, cylindrical needles are recommended; although taper cut or other cutting needles are not contraindicated.
- for cutting the woven graft (ATC) low temperature ophthalmic cautery ($\leq 704^{\circ}\text{C}/1300^{\circ}\text{F}$) is recommended to avoid fraying. If high temperature cautery ($\geq 1204^{\circ}\text{C}/2200^{\circ}\text{F}$) is used, it is recommended to immerse woven prosthesis in sterile saline solution immediately prior to use, preventing focal burning which may occur during cauterization. Do not allow the prosthesis to dry prior to cutting with cautery.
- should the surgical technique require the graft to be clamped, atraumatic clamps (preferably with a piece of soft, or lined jaws) should be applied with the minimum force necessary.

Potential Complications

- Stroke
- Paraplegia
- Paraparesis
- Myocardial Infarction
- Renal Dysfunction
- Ischemia
- Embolism
- Thrombosis
- Bleeding
- Graft Infection
- Wound Infection
- Aneurysm
- Pneumonia
- Amputation
- Death
- Graft dilation
- Graft disintegration
- Graft stenosis
- Pseudoaneurysm
- Intraluminal graft thrombus
- Perigraft air
- Respiratory failure
- Atrial fibrillation
- Chylothorax
- Temporary psychotic syndrome (including transient delirium)
- Hemiparesis
- Visceral ischemia

Precautions and Warnings

- Do not use a prosthesis if the container and/or seal has been opened or damaged, or if the period of

* Not available for sale in the United States or Canada.

sterility has expired.

- The collagen-impregnated graft must never be resterilized.
- Grafts contaminated with blood during the preceding procedures must not be re-used or resterilized.
- The vascular grafts must be handled so as to avoid contact with extraneous particles which, if they adhere to the graft wall, may generate emboli or undesirable interactions with the blood. Furthermore, surgical gloves used to handle grafts should not contain powders, preservatives or lubricants.
- Avoid overstretching the graft; gently expand the graft to smooth the folds.
- Avoid damaging the graft when handling, useatraumatic clamps and appropriate instruments (e.g. vascular clamps). Avoid using these instruments with undue force, otherwise the collagen coating or fabric will be damaged.
- Atraumatic needles are recommended.
- Low temperature ophthalmic cautery ($\leq 704^{\circ}\text{C}/1300^{\circ}\text{F}$) is recommended for cutting woven grafts to avoid fraying.
- AlboGraft Removable Spiral Reinforcement Prosthesis (ASC Models)*: Avoid clamping the graft on its reinforced area.
- AlboGraft Removable Spiral Reinforcement Prosthesis (ASC Models)*: Gently remove the support spiral, otherwise the collagen film will be damaged.
- Care should be taken to ligate and/or cauterize lymphatics in the groin to minimize the occurrence of seroma formation and lymphatic collection subsequent to aorto-femoral or femoropopliteal reconstruction.
- These prostheses should not be implanted in patients who exhibit sensitivity to polyester or materials of bovine origin.

Re-sterilization/Rewrapping

This device is single-use only. Do not reuse, reprocess, or re-sterilize. The cleanliness and sterility of the re-processed device cannot be assured. Reuse of the device may lead to cross contamination, infection, or patient death. The performance characteristics of the device may be compromised due to reprocessing or re-sterilization since the device was only designed and tested for single use. The shelf life of the device is based on single use only. If for any reason this device must be returned to LeMaitre Vascular, place it in its original packaging and return it to the address listed on the box.

Conservation and Storage

The grafts must be kept in their original packaging until used. They must be kept in a dustfree and dry environment and at a temperature between $0^{\circ}\text{C}/32^{\circ}\text{F}$ and $30^{\circ}\text{C}/86^{\circ}\text{F}$.

Safe Handling and Disposal

If serious medical incidents should arise during use of this medical device, users should notify both LeMaitre Vascular and the Competent Authority of the country where the user is located.

This product contains no sharps, heavy metals or radioisotopes, and is not infectious or pathogenic. No special requirements for disposal are evident. Please consult local regulations to verify proper disposal.

Packaging and Shipping of Explanted AlboGraft:

Return of the shipment to LeMaitre Vascular depends on 3 crucial questions:

1. Is the explant obtained from a patient with a known or presumed pathogenic condition at the time of the explant?
2. Is the explant obtained from a patient that has a known treatment history that involves therapeutic radionuclides within the last 6 months?
3. Has the clinician obtained consent from the patient for the sample to be returned to the manufacturer for research purposes?

* Not available for sale in the United States or Canada.

In the event that the response to question 1 or 2 is affirmative, LeMaitre Vascular does not provide adequate guidance for shipment. THESE EXPLANTS SHOULD NOT BE RETURNED TO LEMAITRE VASCULAR UNDER ANY CIRCUMSTANCES. In these cases, the explant should be disposed according to local regulations.

For explants that do not have pathogenic or radiologic hazards, please use the following:

Pre-explantation:

1. If possible, perform a CT or Ultrasound scan of the device to document patency.
2. LeMaitre Vascular can accept clinical information that is patient-anonymized. LeMaitre Vascular requests information including:
 - a) The original diagnosis which resulted in use of the implant,
 - b) The patient's medical history relevant to the implant, including the hospital or clinic at which the device was implanted.
 - c) The patient's implant experience prior to implant removal.
 - d) The hospital or clinic at which the explantation was done, and date of retrieval.

Explantation:

1. Explanted AlboGraft grafts should be directly transferred to a sealable container filled with a solution of alkaline buffered 2% glutaraldehyde or 4% formaldehyde prior to shipment.
2. Cleaning of explanted grafts should be minimal if necessary. Proteolytic digestion should not be used in any circumstances.
3. AlboGraft explants should not be decontaminated under any circumstances. DO NOT autoclave the sample or use ethylene oxide gas to decontaminate.

Packaging:

1. Explants should be sealed and packed in a manner that minimizes potential for breakage, contamination of the environment or exposure to those handling such packages during transit. Material that is absorbent and is cushioning should be selected for isolating the sealable container inside the secondary packaging. Primary and secondary packaging must then be packaged inside an outer package.
2. Explants in sealed primary containers should be labeled with an ISO 7000-0659 Biohazard symbol. The same symbol should be attached to the secondary packaging and to the outer packaging. Outer packaging should also be labelled with Name, Address and Telephone Number of Sender, and the statement, "Upon discovery of damage or leakage, the package should be isolated and sender notified".
3. Packages prepared in the above manner may be shipped to:
LeMaitre Vascular
63 Second Avenue
Burlington, MA 01803

Notices: Limited Product Warranty; Limitation of Remedies

LeMaitre Vascular, Inc., warrants that reasonable care has been used in the manufacture of this device and that this device is suitable for the indication(s) expressly specified in these instructions for use. Except as explicitly provided herein, LEMAITRE VASCULAR (AS USED IN THIS SECTION, SUCH TERM INCLUDES LEMAITRE VASCULAR, INC., ITS AFFILIATES, AND THEIR RESPECTIVE EMPLOYEES, OFFICERS, DIRECTORS, MANAGERS, AND AGENTS) MAKES NO EXPRESS OR IMPLIED WARRANTIES WITH RESPECT TO THIS DEVICE, WHETHER ARISING BY OPERATION OF LAW OR OTHERWISE (INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE) AND HEREBY DISCLAIMS THE SAME. This limited warranty does not apply to the extent of any abuse or misuse of, or failure to properly store, this device by the purchaser or any third party. The sole remedy for a breach of this limited warranty shall be replacement of, or refund of the

* Not available for sale in the United States or Canada.

purchase price for, this device (at LeMaitre Vascular's sole option) following the purchaser's return of the device to LeMaitre Vascular. This warranty shall terminate on the expiration date for this device.

IN NO EVENT SHALL LEMAITRE VASCULAR BE LIABLE FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, SPECIAL, PUNITIVE, OR EXEMPLARY DAMAGES. IN NO EVENT WILL THE AGGREGATE LIABILITY OF LEMAITRE VASCULAR WITH RESPECT TO THIS DEVICE, HOWEVER ARISING, UNDER ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, TORT, STRICT LIABILITY OR OTHERWISE, EXCEED ONE THOUSAND DOLLARS (US\$1,000), REGARDLESS OF WHETHER LEMAITRE VASCULAR HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH LOSS, AND NOTWITHSTANDING THE FAILURE OF THE ESSENTIAL PURPOSE OF ANY REMEDY. THESE LIMITATIONS APPLY TO ANY THIRD-PARTY CLAIMS.

A revision or issue date for these instructions is included on the back page of these Instructions for Use for the user's information. If twenty-four (24) months has elapsed between this date and product use, the user should contact LeMaitre Vascular to see if additional product information is available.

* Not available for sale in the United States or Canada.

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LeMaitre and AlboGraft are registered trademarks of LeMaitre Vascular, Inc.
©2020 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA.
LeMaitre Vascular GmbH